



OSSEointegrated Fixtures

Product Catalogue





Southern Implants is a leading provider of unique and innovative dental implant products with a focus on top-end professional users who want more choices. Southern's expertise in research, development and manufacturing of dental implants allows us to provide Innovative Treatment Solutions that will reduce treatment times and improve patient outcomes.

Striving for excellence and meeting customer needs, has led to our wide product range characterized by Unique and Innovative products which include;

- Multiple interfaces, to suit customer preference.
- INVERTA™ implant, featuring a body-shift design, engineered for primary stability and suitable for immediate loading.
- Co-Axis®, sub-crestal angle correcting implant, available in angulations of 12, 24 & 36° and various internal and external connections.
- MAX implant, specifically designed for immediate molar tooth replacement.
- The ZYGAN and ZYGELEX implants for severely resorbed maxilla and craniofacial reconstruction.

Our product portfolio is in synchronized evolution with protocol improvements and technological advances.

My sincere thanks to all specialists, dentists and technicians who put their trust in our company.



Graham Blackbeard
Managing Director, Southern Implants

CONTENTS

External Hex Ultra-Short Ø3.75mm / Ø4.5mm Implant

Implants & Surgical Components.....	Page 06
Prosthetic Flowchart.....	Page 07
Prosthetic Flowchart: Magnetic Retention, Steco Titanmagnetics® for IE.....	Page 07
Site Preparation Sequence (Parallel walled & Tapered implants).....	Page 08
Instrumentation (Parallel walled & Tapered implants).....	Page 09
Surgical Tray	
I-EO-EG.....	Page 10
Driver and Screw Information	
.....	Page 12
Steco Titanmagnetics	
.....	Page 12
Implant dimensions and Information	
.....	Page 14
Explanation of Symbols	
.....	Page 15

NOTE:


- Images are for illustration purposes only and do not necessarily accurately represent the product.
- All dimensions in this catalogue are in mm, unless otherwise specified.
- Not all products are cleared for sale in all countries.

EXTERNAL HEX

**Ø3.75mm
Ø4.5mm**

ULTRA-SHORT

Ø3.75mm Implants (Parallel Walled)

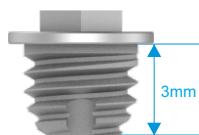**IE3**

Implants are available in lengths of:

Implant dimensions and information - page 14

Cylindrical

Item code	Implant lengths (in mm)
IE	3

Surgical Components





Cover Screw



Healing Abutments

TBE8n

Ø3.75mm / Ø4.5mm Implants (Parallel Walled & Tapered)

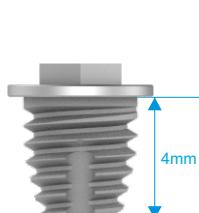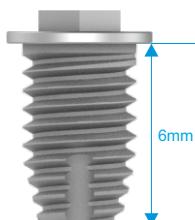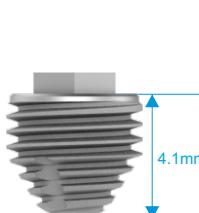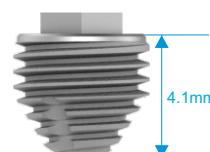**IE4****IE6****IET4****IETi4**

(Internal Drive)

Implants are available in lengths of:

Implant dimensions and information - page 14

Cylindrical

Item code	Implant lengths (in mm)
IE	4 / 6

Tapered

Item code	Implant lengths (in mm)
IET	4
IETi	4 (Internal Drive)

Surgical Components


Cover Screw



Healing Abutments

SCU22/3/4/5/6/8
lengths***TB**2/3/4/5/6/8
lengths***WB**2/3/4/5/6/8
lengths

Prosthetic Flowchart

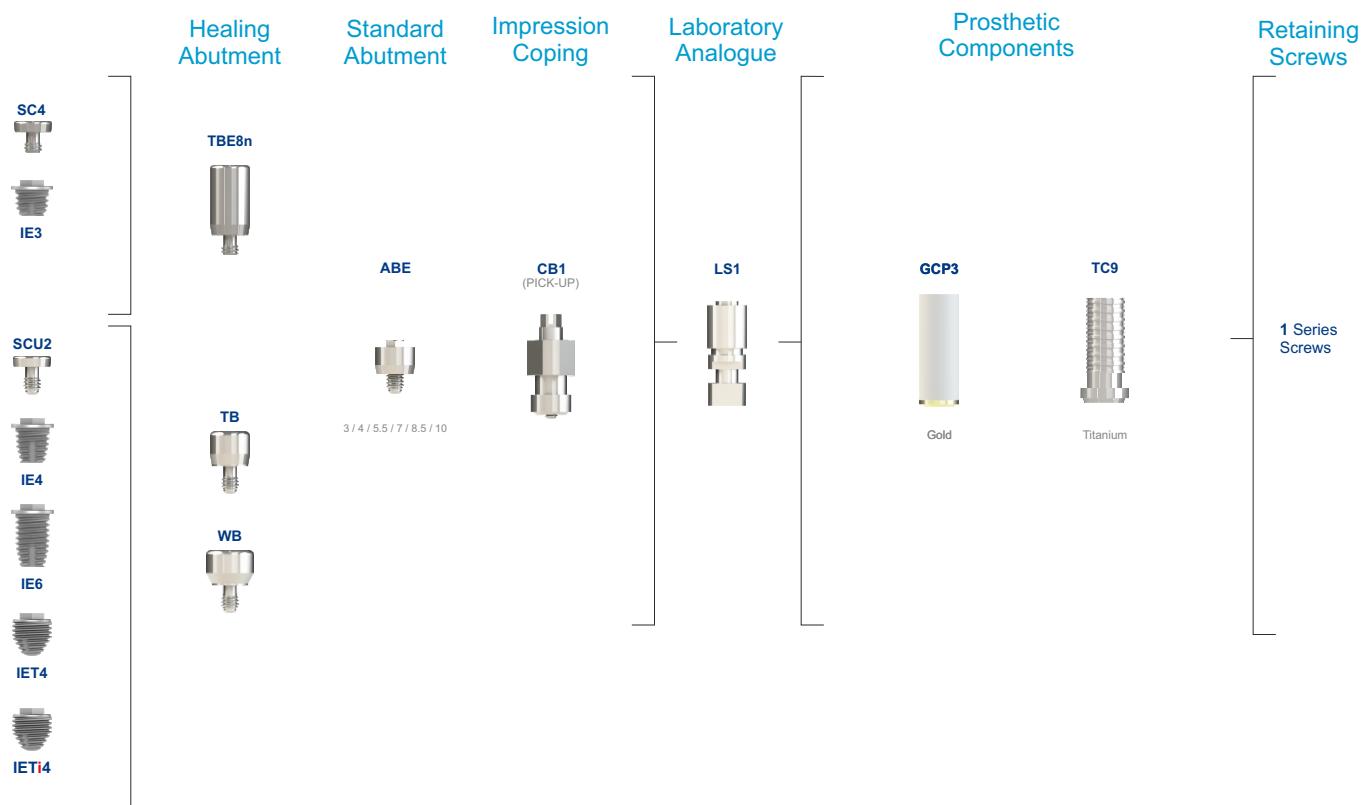

Prosthetic Flowchart

Magnetic Retention, Steco Titanmagnetics® for IE




NOTE: See page 20-23 for more information on Steco Titanmagnetics® prosthetics and instrumentation.

Site Preparation Sequence

Ø3.75mm Parallel Walled (IE3)

Soft Bone


D-RB-MS
Ø2.0mm

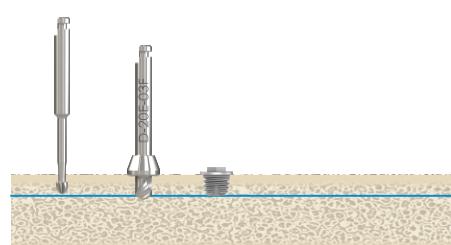

Medium Bone

D-RB-MS
Ø2.0mm

Ø3.0mm

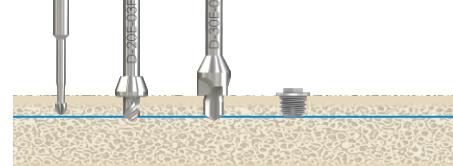


Dense Bone

D-RB-MS
Ø2.0mm

Ø3.0mm

D-TAP-BS

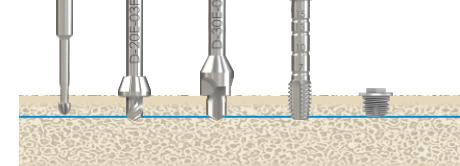


3mm length

Ø3.75mm Parallel Walled (IE4)

Soft Bone


D-RB-MS
Ø2.0mm

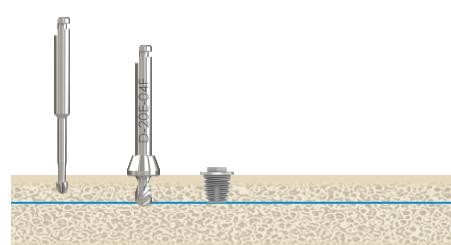

Medium Bone

D-RB-MS
Ø2.0mm

Ø3.0mm

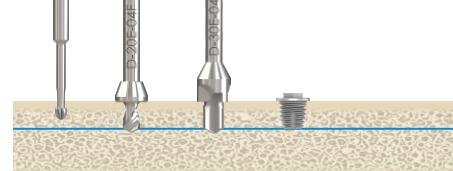


Dense Bone

D-RB-MS
Ø2.0mm

Ø3.0mm

D-TAP-BS

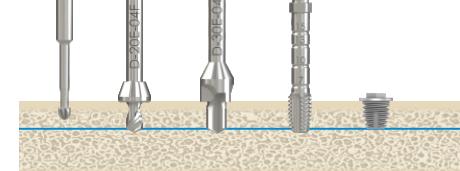


4mm length

Ø3.75mm Parallel Walled (IE6)

Soft Bone


D-RB-MS
Ø2.0mm

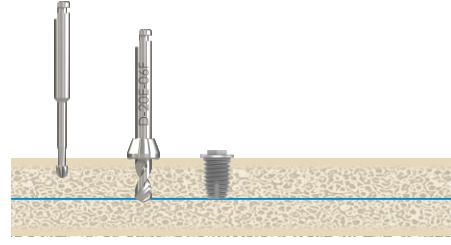

Medium Bone

D-RB-MS
Ø2.0mm

Ø3.0mm

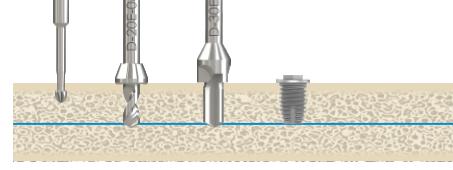


Dense Bone

D-RB-MS
Ø2.0mm

Ø3.0mm

D-TAP-BS

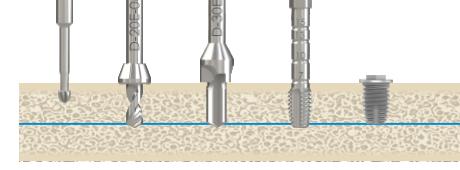


6mm length

Ø4.0mm Tapered (IET4 / IETi4)

Soft Bone


D-RB-MS
Ø2.0mm

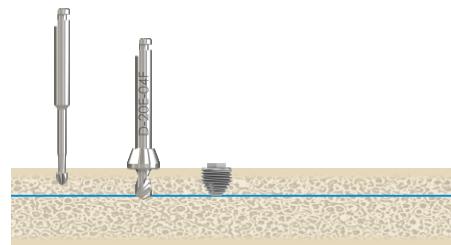

Medium Bone

D-RB-MS
Ø2.0mm


Ø4.0mm

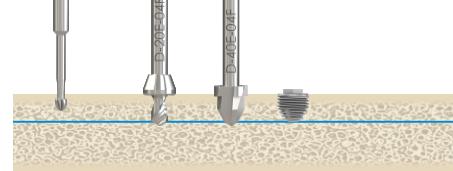

Dense Bone

D-RB-MS
Ø2.0mm

Ø4.0mm

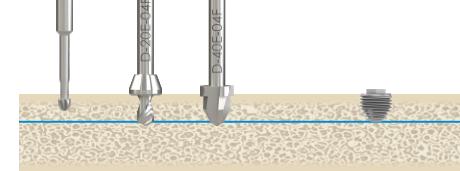

4mm length

NOTE: The total implant length is reached when the flange of the 4mm tapered drill touches the bone surface. Drilling beyond this point may compromise primary stability of the implant.

INSTRUMENTATION

Fixture Mounted Implants (Parallel Walled & Tapered)

	Insertion Tools				Wrench Converters	
	I-CON-X	I-CONU-X	I-WI-CST	I-WI-SH		
						
IE3	Fixture Mount Driver (With W&H)	Fixture Mount Universal Driver (Without W&H)	For Handpiece inserts (Latch-type) featuring the W&H hex.	For HEX connection of fixture mounts		
						
IE4						
						
IE6						
						
IET4						


Internal Drive Implant (Tapered)


Ø4.5mm Tapered (**IETi4**)

	Insertion Tools			Wrench Converters	
	I-HID-S / M / L	I-HIDU-S / M / L	I-WIID-S / L	I-WI-CST	I-WI-SS
					
IETi4 (Internal Drive)	Handpiece Insert with W&H	Handpiece Insert Universal	Wrench Insert	For Handpiece inserts (Latch-type) featuring the W&H hex.	For SQUARE connection of fixture mounts and instruments

* Most instruments are available in various lengths

I-EO-EG For surgical placement of Extraoral Implants

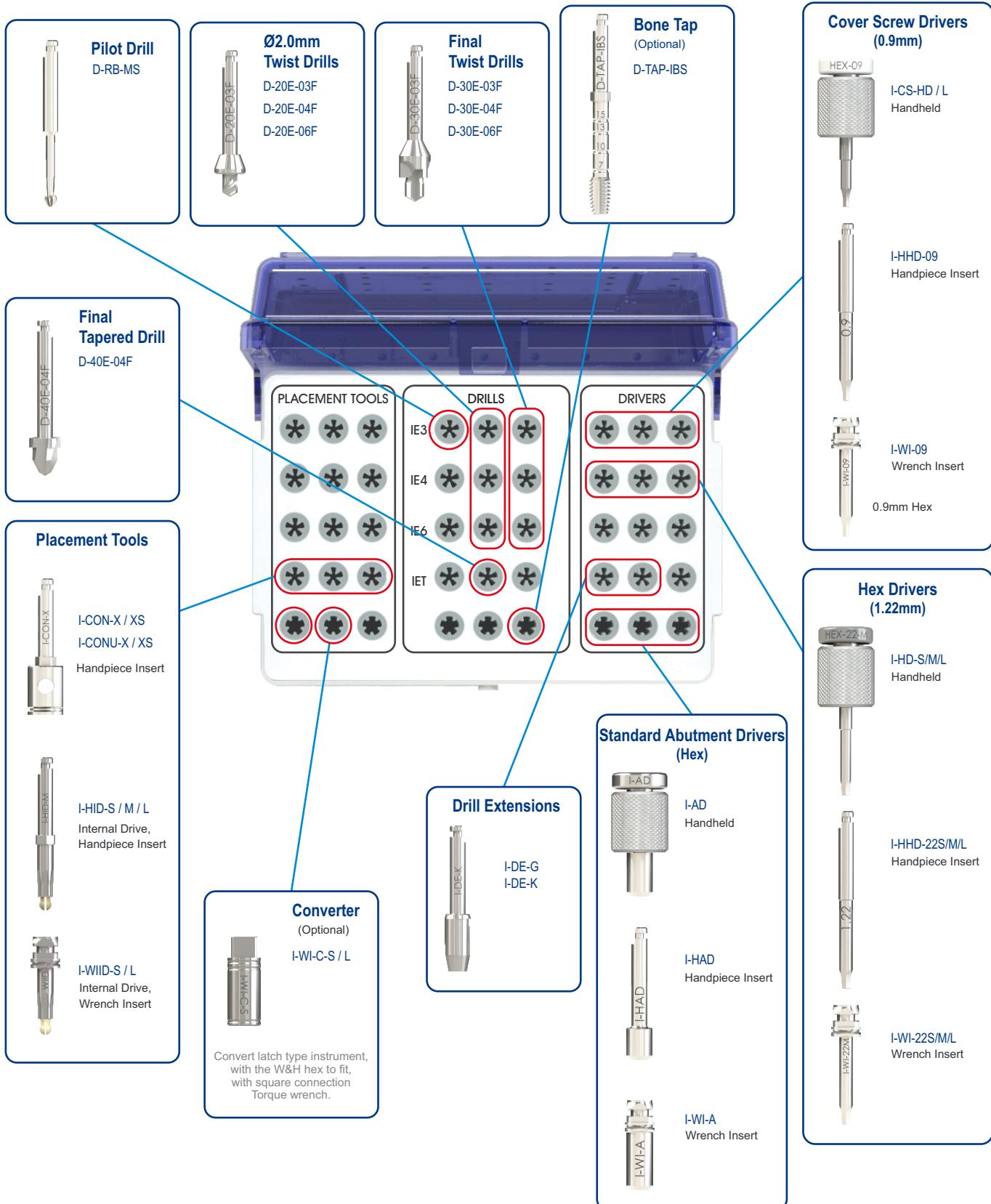

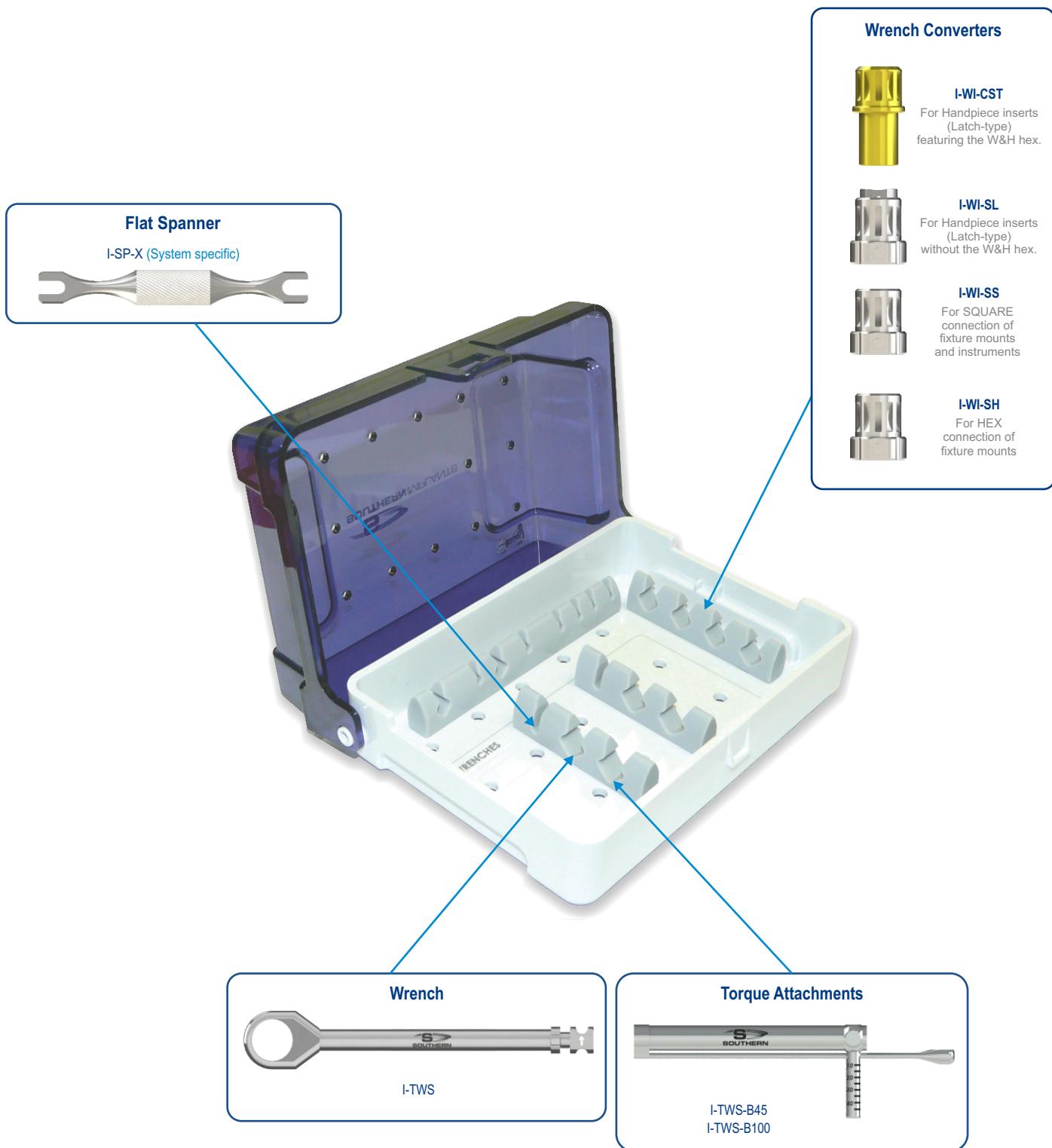

NOTE:

- The surgical kit has an intuitive layout to guide the surgeon through the drill sequence.
- Most Instruments available in Short / Medium / Long.
- All instruments and tooling used during the procedure must be maintained in good condition, and cleaned and sterilized prior to use. Please consult the Southern Implants Cleaning and Sterilization Procedure Guidelines (CAT-1039) for guidance concerning the maintenance of drills, instruments, and surgical trays.
- For Polishing Protector Caps refer to CAT-1010

DRIVER AND SCREW INFORMATION

Alternative Drivers

Hex Drivers



NOTE: Refer to CAT-1203 for alternative drivers

1 Series Screws



Screw Head Connections

Laboratory Screw

1.22 Hex

BSH1*
(Blackened and for laboratory use only)



Hex



Unigrip


NOTE:

- Due to design revisions screw tips may be flat or rounded.
- Always ensure that the correct screw is used for the relevant implant and component.
- Refer to CAT-1003 for alternative slotted 1 Series screws

STECO TITANMAGNETICS


Magnetic attachments for facial prosthesis

The universal and easily manageable solution for implant-supported face prostheses and obturators.

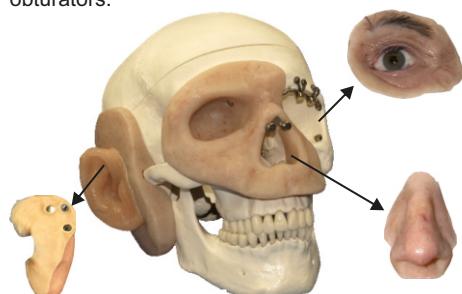


Dual magnet system

Two open magnet fields exert a force of attraction over a longer distance than mono-magnet systems.

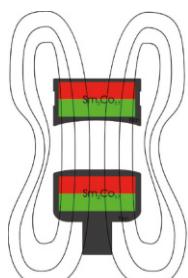


The construction principle

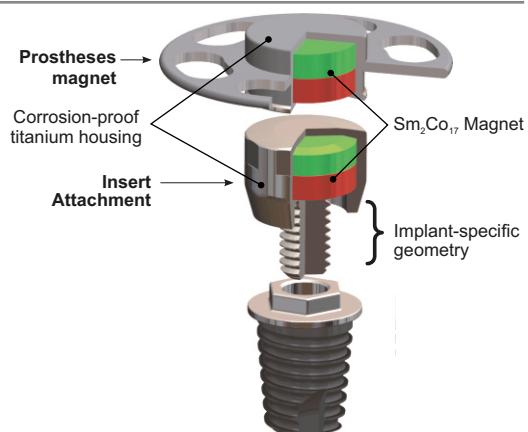

Telescopic magnets


The StecoTitanmagnetics® **T-Line** was developed exclusively for extraoral use. The telescopic bearing of 2.5 mm gives the facial prosthesis secure fixation on the tissue. The design allows a certain amount of rotation and axial movement, without losing the withdrawal force.

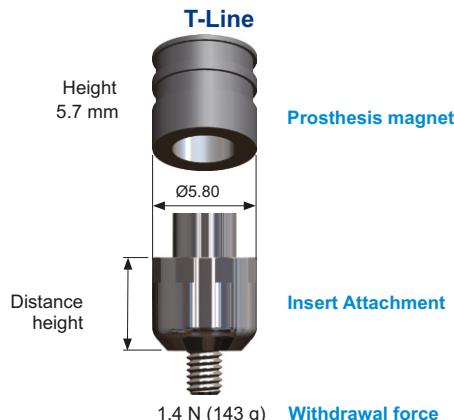


Spherical magnets

The slightly curved surface of the spherical StecoTitanmagnetics® X-Line and the Z-Line facilitates self centering and has the lowest lateral force introduction into the implant. X-Line and Z-Line are therefore particularly suitable for short implants. Due to the absence of lateral steering, they are also indicated for highly divergent implants.

X-Line

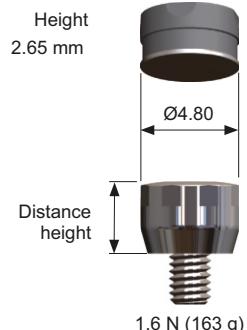

Z-Line

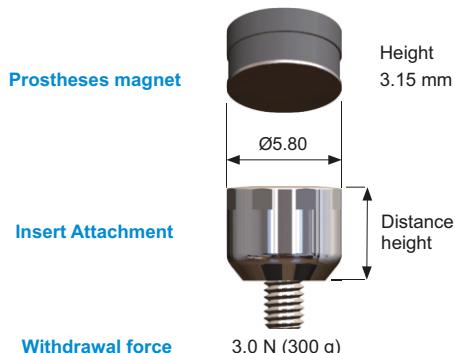

Prosthetic Components

In X-Line and Z-Line there are prostheses magnets for acrylic and, with additional retention, for silicone. For more lateral support there are prostheses magnets with extended collars available.

For silicone



For acrylic



In the T-Line there are also prostheses magnets for acrylic and with additional retention for silicone.

For silicone



For acrylic



Accessories

To make a new prostheses, each product line has its own compatible impression copings and laboratory analogues.

Impression Copings

- magnetic impression without screws
- lock magnetically to the insert
- the external geometry ensures the impression material is securely fixed

Laboratory Analogues

- for making models quickly and hygienically
- the original implant does not have to be used on the insertion point of an analogue
- there is no need for time-consuming cleaning
- a length of 16 mm is also available on request



Healing Flange

The Healing-Flange is placed on the insert to stabilize the bandage on fresh skin perforation tissue, after application of the insert. This minimizes excessive tissue proliferation. The Healing-Flange has its own magnet core and is made off tissue friendly Titanium.



Insertion standard tools

X-Line


H.00.04.X2	
Connector to handpiece	
H.00.04.X4	
Connector to ratchet	
H.00.04.X1	
Connector to I-Ratchet-2	

Z-Line & T-Line

H.00.04.T2	
Connector to handpiece	
H.00.04.T4	
Connector to ratchet	
H.00.04.T1	
Connector to I-Ratchet-2	

IMPLANT DIMENSIONS AND INFORMATION

RANGE										IMPLANT LENGTHS		
	MAJOR DIAMETER	PLATFORM DIAMETER (flange)	Prosthetic Colour Code	HEX WIDTH	HEX HEIGHT	COLLAR HEIGHT	THREAD PITCH	APEX DIAMETER	Cylindrical or Tapered C/T	3	4	6
IE	Ø3.75mm	3.75	4.80	2.70	0.7	0.40	0.6	2.80	C	✓	✓	✓
IET	Ø4.5mm	4.5	4.07	2.70	0.7	0.25	0.5	2.45	T		✓	
IETi	Ø4.5mm	4.5	4.07	2.70	0.7	0.25	0.5	2.45	T		✓	

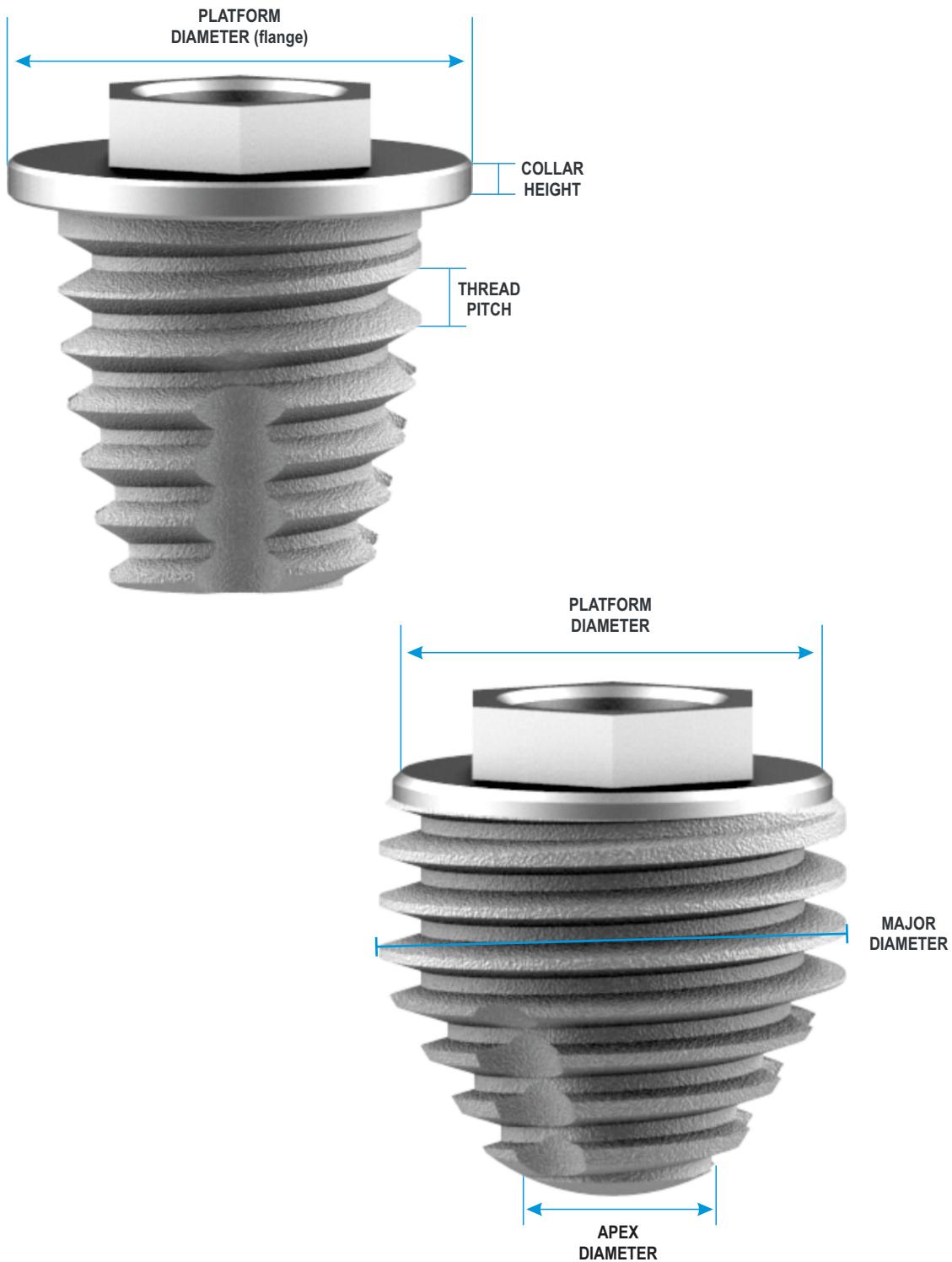

EXPLANATION OF SYMBOLS

The following symbols are used on our packaging labels and they indicate the following:

- 1  Manufacturer
- 2  Implant image
- 3  Implant details and size
- 4  Sterilization using Irradiation



Do not Resterilize



Consult instruction for use



Do not reuse



Caution



CE mark and notified body number



Expiry date





Sterile unless package is opened or damaged

- 5  2D Bar coding
Contains the GTIN, Expiry Date and LOT Number
- 6  Patient sticker for documentation purposes
(to be used by health care provider on patient file)

- 7  Prescription device

**CAUTION: FEDERAL LAW RESTRICTS THE DEVICE TO SALE BY
OR ON THE ORDER OF A LICENCED HEALTH CARE PROVIDER.**

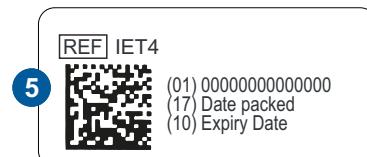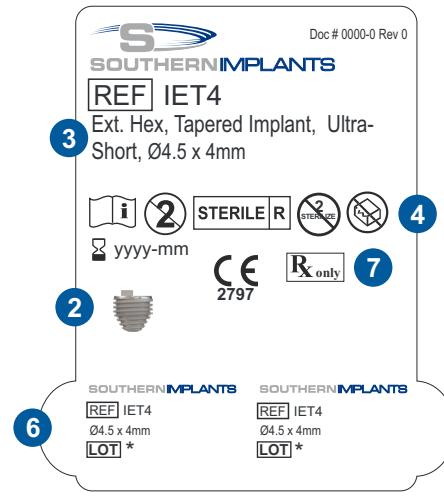

For more information, please contact your
Southern Implants Representative or visit southernimplants.com



Subsidiaries

Australia

Southern Implants Australia
T: +61-2-8076-9337
E: info@southernimplants.com.au

United Kingdom

Southern Implants UK
T: +44-20-8899-6845 / 6 / 7
E: info@southernimplants.co.uk

USA and Canada

Southern Implants North America Inc.
T: +1-561-472-0990
E: customercare@southernimplants.com

Southern Implants Products are distributed world-wide, please visit southernimplants.com for a list of Distributors.

©2019 Southern Implants (Pty) Ltd. All rights reserved. Co-Axis, Southern Implants, and Subcrestal Angle Correction are registered trademarks and Body-Shift, INVERTA and PROVATA are trademarks of Southern Implants (Pty) Ltd. Not all implants are available in all markets.